

USER GUIDE

User Guide for the Laue-Camera Application

SOFTWARE-VERSION: 1.1.1
DATUM: 19.01.2022
AUTHOR: JULIAN KAISER

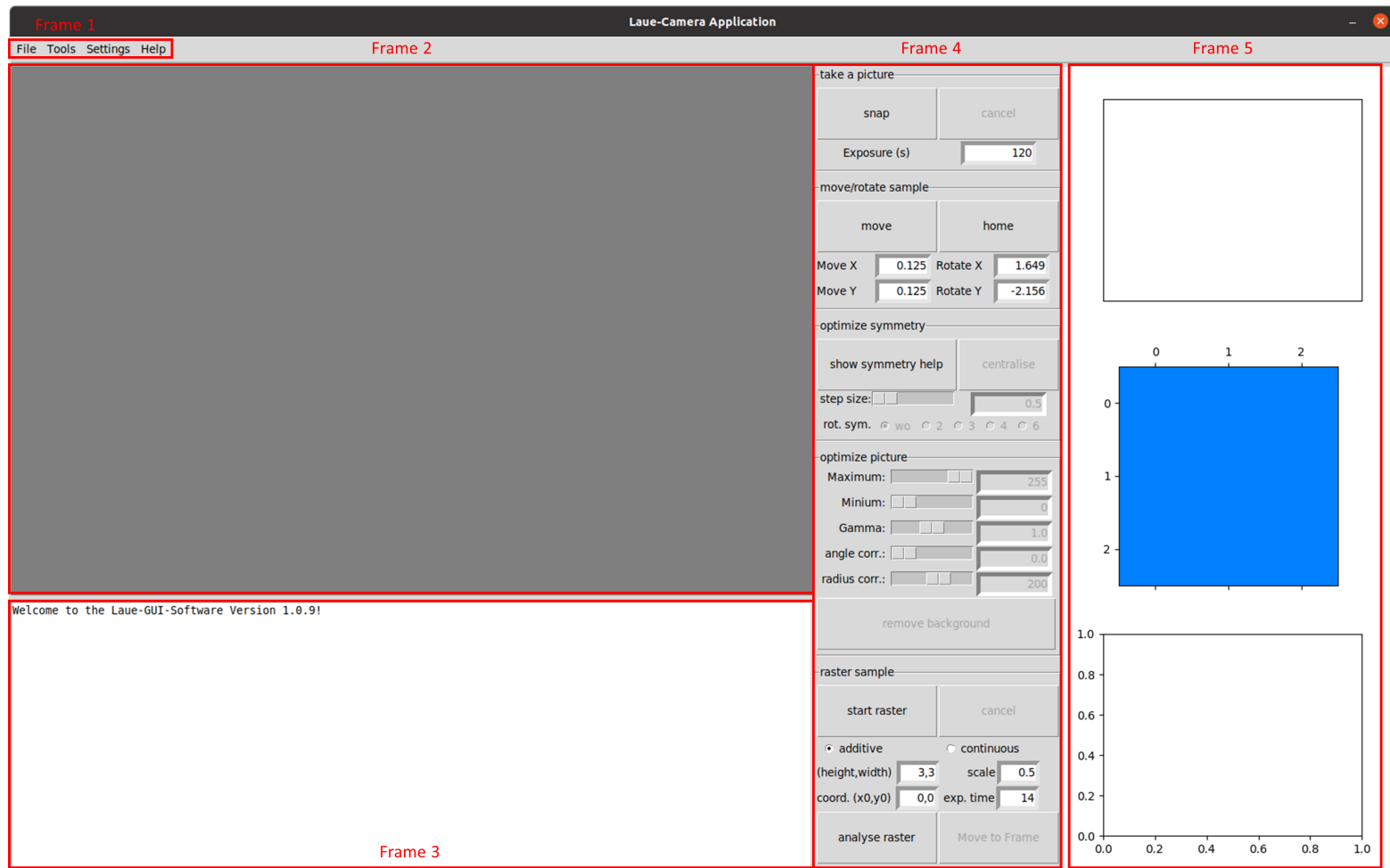


Fig. 0.1: Laue-Camera Application after starting

Contents

System requirement	4
Changelog	5
1.1.1	5
Frame 1 - menu	6
File	6
Tools	7
Settings	12
Frame 2 - image viewer	17
Frame 3 - console	17
Frame 4 - area of operations	19
Section 1 - image acquisition	19
Section 2 - translation and rotation of the sample	19
Section 3 - symmetrization	20
Section 4 - optimize picture	22
Section 5 - raster sample	28
Frame 5 - diagram viewer	30
Diagram 1 - maximum finder	30
Diagram 2 - contour diagram	30
Diagram 3 - backgroundfit	31
Outlook	32

System requirement

This application is coded in python version 3.8.10 for PSLViewer version 3I-CFL from Photonic Science. It's purpose is to simplify and optimize the handling of the Laue-Diffraktometer from the chair of applied physics 6, the motorized goniometer and the image processing. Basic functions for the communication with the server and goniometer are adopted from Marvin Klinger and enhanced.

It is recommend to use the Laue-Diffraktometer and this application only with 20 kV and 30 mA, because this application was developed specifically for this values. Furthermore it is optimized for the visualisation of Na_2IrO_3 crystals, which produce rather bad images. For other crystals, all the possible functions for image processing may don't be necessary. Also other settings may be more beneficial, so many default settings can be changed. If any non changeable value should be changeable, please contact the author! Moreover not only the goniometer is accessible, but also the Windowscomputer with the corresponding programm PSLViewer, so that a full remote control is possible. If PSLViewer crashed, two important settings on the Windowscomputer have to be adjusted in advance, that a remote restart is possible:

- The application „cports.exe“ has to be in the folder with the same name at „Local Disk (C:)“. If it was deleted, a download is possible [here](#).
- At „Task Scheduler Library“ the task „restart_PSL_Viewer_2“ has to exist. If not, this [tutorial](#) will explain the required steps (without changing the icon). The application that has to be selected is „restart_PSL_Viewer_2.bat“ at the folder „PSLViewer-3I“ and the resulting connection must be callable from „restart_PSL_Viewer.bat“ at the folder „PSLViewer-3I“.

Changelog

1.1.1

- Convert a .tif- to a .bmp-raster
- Detection of areas with complete different alignments through the detection of the center and the orientation of the diffraction pattern
- Different visualization options for rasters with center detection data
- New evaluation of a kontur for long exposure times
- Save a flipped image at the horizontal or vertical image center

Frame 1 - menu

At the menu bar in the sub menu [File](#), the user can save the displayed picture on [Frame 2](#) and display a recorded image or raster. In the sub menu [Tools](#), several tools can be used and in [Settings](#) the user can change the value of important parameters. In the sub menu [Help](#), this PDF file will be displayed.

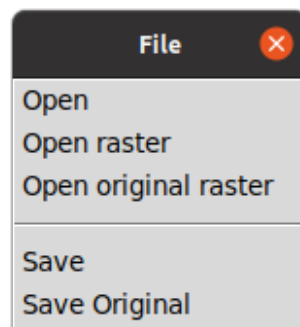


Fig. 1.1: File-menu with all available functions

File

- **Open**: Opens a dialog to select a recorded grayscale image with the dimensions 975x643 pixel, either a 8-bit („.bmp“) or higher data formatted („.tif“) images. This command is also available with the shortcut „STRG + o“.
- **Open raster**: Opens a dialog to select a recorded raster, saved as 8-bit and multiplot image. The dialog only displays folders, so the user has to select only the corresponding folder (dialog has to be inside the folder). For correct visualization, the order of the saved rasterframes is very important (name-sorted must be equal to the acquisition order). This command is also available with the shortcut „STRG + r“.
- **Open original raster**: Opens a dialog similar to „**Open raster**“ with the same requirements on the file-order, but in this command only „.tif“ are opened. This command is also available with the shortcut „STRG + SHIFT + R“.
- **Save**: Opens a dialog to save the displayed 8-bit image at [Frame 2](#) as „.bmp“ file. This command is also available with the shortcut „STRG + ALT + s“. Additionally, the shortcut „STRG + s“ also saves the image, but doesn't open a dialog but

saves the image at the selected directory in [Settings](#). The file name contains the name/number of the sample, applied voltage and current, distance to the detector, exposure time, comment and a uprising number starting at „000“, all seperated by „-“ (more information to the parameters is available at [save settings](#)).

- **Save Original:** Opens similar to „Save“ a dialog to save the original image (greater than 8-bit). This command is also available with the shortcut „STRG + ALT + SHIFT + S“. Also the shortcut „STRG + SHIFT + S“ saves similar to „Save“ the original image at the selected directory in [Settings](#) with the suffix „original“ in the filename.

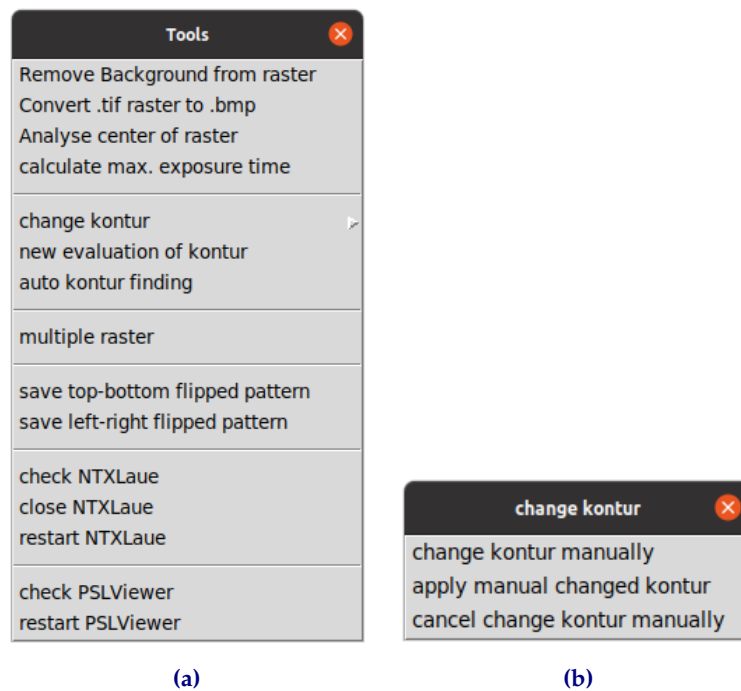


Fig. 1.2: Tools-menu with all available functions

Tools

- **Remove Background from raster:** Opens a dialog similar to „[Open raster](#)“ to select a original raster (no multiplot images!) to remove the background of each frame. After the selection, a new dialog opens to select the save directory, where a new folder is created for the new calculated raster. The command subtracts the background with the selected options at [Settings](#) from the whole image.
- **convert .tif to .bmp raster:** Converts a original raster (.tif files) to a raster with 8-bit (.bmp files). Therefore the user has to set a mean value and then the maximum

pixel value is tuned, so that the mean value of the resulting image is near to the defined one. The mean value can be set in a range from 255 (maximum) to the mean value of the original image (minimum).

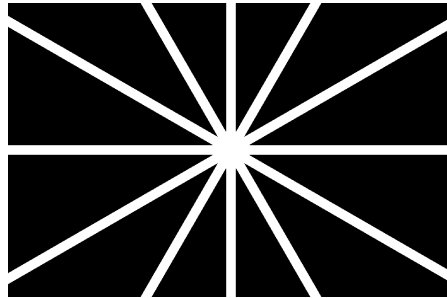


Fig. 1.3: Example of a mask with six lines and a linewidth of 21 pixel

- **Analyse center of raster:** Opens a dialog to select a raster to analyse and afterwards a dialog to select the save directory for the collected data. The method tries to detect the center and the smallest clockwise degree a „line“ is rotated from the vertical image-center of the centered diffraction image. Therefor, according to the [Settings](#), a mask slightly larger as the image is generated first (see fig. 1.3). The center of this mask is put on each pixel of the search area and then the mean value of all pixel who aren't covered by the mask is calculated. To define the rotated degree, the mask is additionally rotated in steps of 1° until it matches with the original one and for each orientation the mean value of the remaining pixel is calculated. The idea behind the method is, that for the best orientation of the mask on the image the mean value of the remaining pixel is at a minimum. First for each angle of rotation the mean value of the mean values from the whole search area is calculated and the smallest one is interpreted as the best orientation. Second the center of the diffraction pattern is the center of mass from all minima of the mean values in the search area from the corresponding angle. The results are saved in the logfile from the selected raster additionally. Furthermore if a raster with data from the center detection method is loaded, the results can be visualized in [Frame 5](#). With a right-click on the kontur the user can choose between the normal kontur, the different angle of rotations (each color represents one specific angle from minimum (dark blue) to maximum (dark red)), the mean distance from the center to surrounding centers of the diffraction images and the combination of the angle of rotation and distance map (added together).
- **calculate max. exposure time:** Calculates the new maximum exposure time. More information are available at [Section 1](#) in [Frame 4](#)
- **change kontur:** In this sub menu, the user can change an acquired contour, if the algorithm misvalue a frame in a raster. The command „**change kontur manually**“ makes the current contour at [Frame 5](#) changeable and the user can navigate through the contour with the arrow keys. To change the selected frame, the user

has to type in numbers from 0 to 3. With the command „**apply manual changed kontur**“, a dialog opens to select a logfile.txt, which contour will be overwritten with the displayed contour. With the command „**cancel change kontur manually**“ the displayed contour isn't changeable any more.

- **new evaluation of kontur**: Opens a dialog to select a .tif-raster, from which the kontur will be reevaluate. This method is necessary for long exposure times, because for an appropriate characterization the remove of the background is necessary. The remove of the background is done by the selected method in [Settings](#), rescaled to an 8-bit image and then evaluated by the same algorithm as in a normal kontur acquisition. The new results are saved in the logfile of the selected raster.

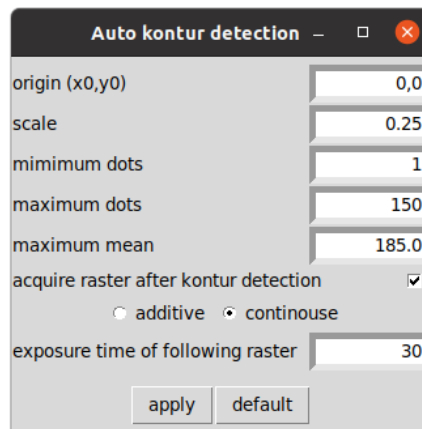


Fig. 1.4: „auto kontur finding“ Dialog with all available functions

- **auto kontur finding**: Opens a dialog to configure the search criteria:
 - **origin** (x_0, y_0): origin of contour search; $x_0, y_0 \in \mathbb{R}$
 - **scale**: step size between two frames; $\text{scale} \in \mathbb{R}^+$
 - **minimum dots**: minimum number of detected points, that are necessary to interpret the frame as a hit; $\text{min}_{\text{dots}} \in \mathbb{N}$
 - **maximum dots**: maximum number of detected points, that are required to interpret the frame as a hit; $\text{max}_{\text{dots}} \in \mathbb{N}$
 - **maximum mean**: maximum mean value up to a frame is interpreted as a hit; $\text{max}_{\text{mean}} \in \mathbb{R}^+$
 - **acquire raster after kontur detection**: acquire after the detection of the contour a new raster with the determined parameters;
 - **additive/continuous**: the user can change between additive (max. 16 s exposure) or continuous (max. exposure) acquisition (for more information see [Section 1 in Frame 4](#))

- **exposure time of following raster**: exposure time of a frame in the subsequent raster, $\text{exposure} \in \mathbb{N}$

This method tries to determine the contour of the sample by moving the sample from the origin in steps of step size in every direction until no diffraction pattern can be detected anymore (number of dots isn't in the defined range). Afterwards all edges are analysed to ensure that the whole sample is inside the contour. The exposure time is 31 s.

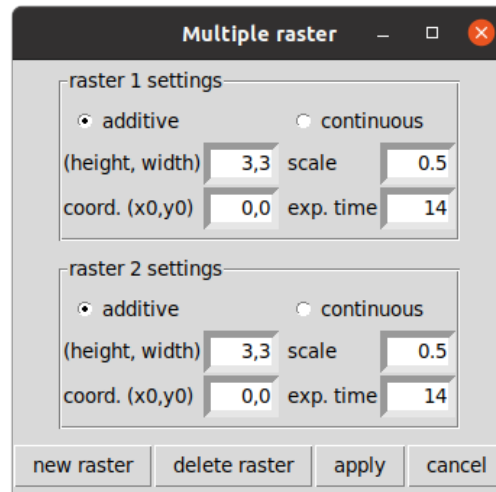


Fig. 1.5: „multiple raster“ Dialog with all available functions

- **multiple raster**: Allow the user to acquire from the same sample various raster. To configure all settings, a dialog will open where more than one raster can be selected with the same settings as a normal raster in [Section 5](#). With the button „apply“, all selected rasters are acquired in a row.
- **save top-bottom flipped pattern**: Saves the image from the diffraction pattern mirrored to the horizontal center of the image because of the same reason in **save left-right flipped pattern**. Is the sample inserted with left and right match with the left and right from the detector, so top and bottom is inverted due to the 180° rotation. Mirror on the horizontal or vertical image center depends therefor on the insertion of the sample inside the goniometer.
- **save left-right flipped pattern**: Saves the diffraction pattern mirrored to the vertical center of the image. This can be helpful for the problem of a correct direction determination on a sample, which results of the different point of views on the image from the microscope and the diffraction pattern. The point of view through the microscope (eye) at the sample is mirrored to the point of view from the computer to the detector (see [fig. 1.6](#)). Therefore one image has to be mirrored vertically for this insertion of the sample.

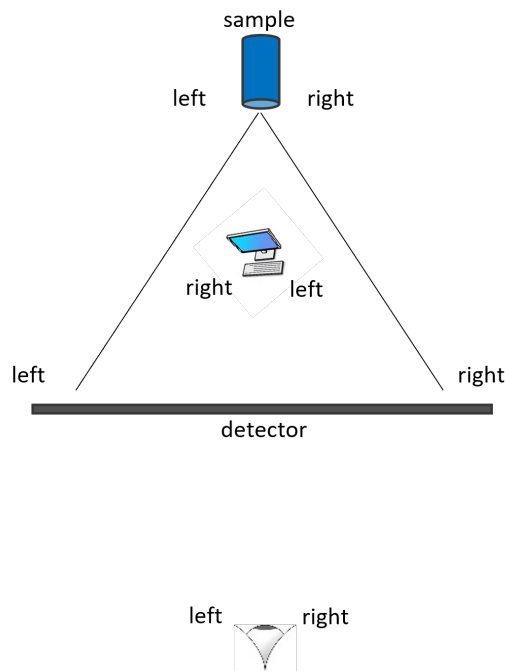


Fig. 1.6: Visualization of the necessity to mirror the diffraction pattern or the image from the microscope for a correct determination of the direction in the sample. Rotating the sample to match right an left lead to swapping top and bottom.

- **check NTXLaue:** Checks the current status of the camera (if the application is opened at the other computer) and displays the status in [Frame 3](#).
- **close NTXLaue:** Closes the camera and displays the status in [Frame 3](#)
- **restart NTXLaue:** Restarts the camera and displays the status in [Frame 3](#). If the camera is closed, the camera will be opened.
- **check PSLViewer:** Checks if the application „PSLViewer“ and the data connection to „server control“ is open and connected.
- **restart PSLViewer:** Restarts the application „PSLViewer“ with „server control“ on the other windows computer with a batch-file. It is useful, if the application crashed or isn't responding. To restart the application, the data connection still has to be open to send the command to the other computer. With this command, a reasonable remote control of the PSLViewer application is possible via, for example, TeamViewer.

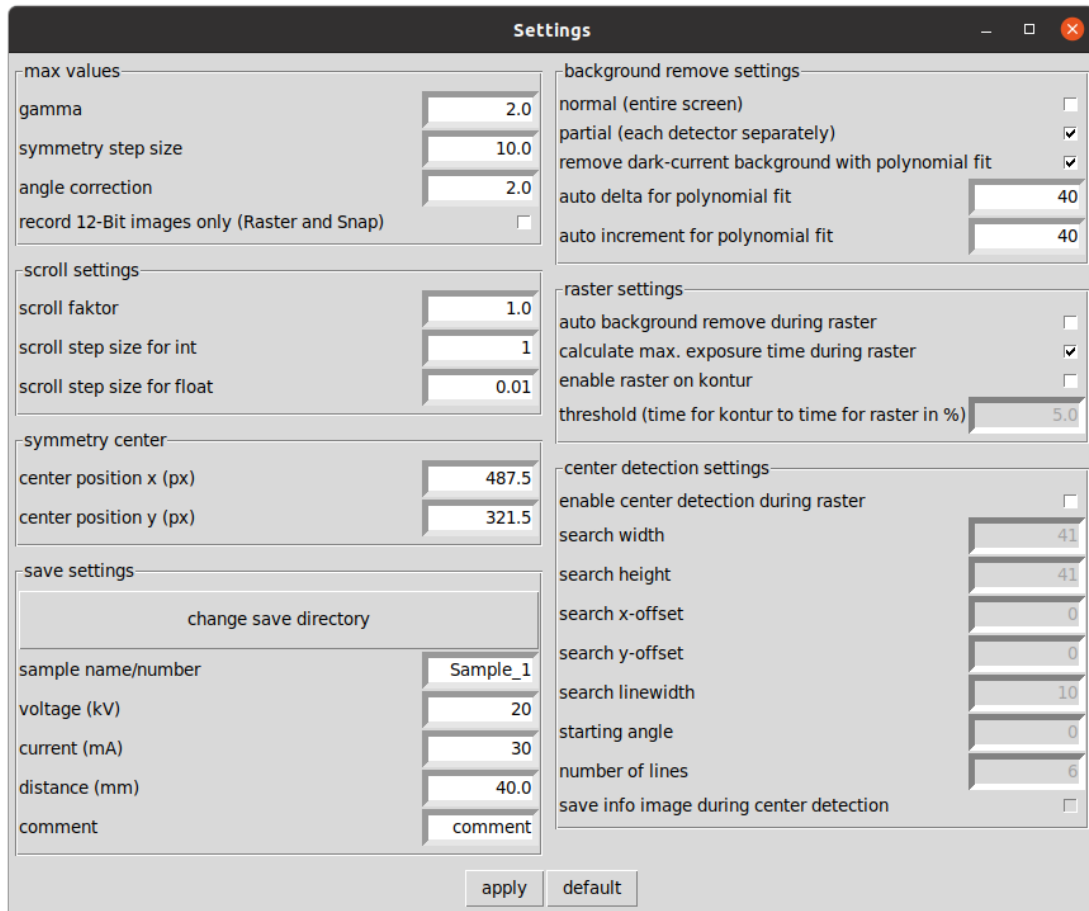


Fig. 1.7: „Settings“ Dialog with all available functions

Settings

max values:

- **gamma**: maximum value of gamma (Minimum is 0); $\max_{\text{gamma}} \in \mathbb{R}^+$
- **symmetry step size**: maximum value of the angle of rotation for the symmetry help (Minimum is 0.5); $\max_{\text{sym}} \in \mathbb{R}^+$
- **angle correction**: maximum value of the angle of rotation for the „angle correction“ (Minimum is 0); $\max_{\text{angleCorr}} \in \mathbb{R}^+$

scroll settings:

- **scroll faktor**: factor which will be multiplied to the stepsize of the sliders in [Frame 4](#); $\text{faktor} \in \mathbb{R}^+$

- **scroll step size for int**: stepsize of the sliders in [Frame 4](#) with integers; $\text{stepsize}_{\text{int}} \in \mathbb{N}$
- **scroll step size for float**: stepsize of the sliders in [Frame 4](#) with floats; $\text{stepsize}_{\text{float}} \in \mathbb{R}^+$

symmetry center:

- **center position x (px)**: X-coordinate of the point, on which the diffraction pattern will be symmetrized in [Frame 4](#) in [Section 3](#). Unfortunately, not all points on the detector are suitable for a center position due to the limitation of the goniometer. Is state in pixels between 0 and 975. $\text{center}_x \in \mathbb{R}^+$
- **center position y (px)**: Y-coordinate of the point, on which the diffraction pattern will be symmetrized in [Frame 4](#) in [Section 3](#). Unfortunately, not all points on the detector are suitable for a center position due to the limitation of the goniometer. Is state in pixels between 0 and 643. $\text{center}_y \in \mathbb{R}^+$

save settings:

- **change save directory**: Opens a dialog to select a directory, which is the origin for save and load dialogs and the save directory for the auto save commands.
- **sample name/number**: Defines the name or number of the current sample, which is used for the filename in the auto save commands. Shouldn't contain the symbol „/“, due to possible errors while saving.
- **voltage (kV)**: Defines the applied voltage at the x-ray tube in kV. Attention: other values than the default value may result in bad or false results during background subtraction! $U \in \mathbb{N}$.
- **current (mA)**: Defines the applied current at the x-ray tube in mA. Attention: other values than the default value may result in bad or false results during background subtraction! $I \in \mathbb{N}$.
- **distance (mm)**: Defines the distance between sample and detector. Value is also used in [Frame 4](#) in [Section 3](#) to calculate the rotational angles to symmetrize the diffraction pattern. $d \in \mathbb{R}^+$
- **comment**: Defines a addition comment which is added to the filename in the auto save commands. Shouldn't contain the symbol „/“, due to possible errors while saving.

background settings:

- **normal (entire screen)**: Fitfunction for background remove is calculated for the whole area of the detector (see also [remove background](#))

- **partial (each detector separately)**: Fitfunction for background remove is calculated and removed for every half of the detector separately (see also **remove background**)
- **remove dark-current background with polynomial fit**: If activated, it removes also the background from the dark-current with a polynomial function (see also **remove background**)
- **auto delta for polynomial fit**: Defines the parameter δ for the automatic calculation of the maximum pixel value in the subtraction of the background due to the dark-current. (see also **remove background**). $\delta \in \mathbb{N}^+$
- **auto increment for polynomial fit**: Defines the parameter inkrement for the automatic calculation of the maximum pixel value in the subtraction of the background due to the dark-current. (see also **remove background**). Inkrement $\in \mathbb{N}^+$

raster settings:

- **auto background remove during raster**: If selected, the background of each frame from every raster with an exposure time greater than 60 s will be subtracted while the raster is acquired. Is saved in an addition folder.
- **calculate max. exposure time during raster**: If selected, the new maximum exposure will be calculated on every frame.
- **enable raster on kontur**: If selected, before a raster is acquired, a faster raster with 31 s exposure time will be acquired and according to the resulting contour, only the frames where a hit was detected will be acquired.
- **threshold (time for kontur to time for raster)**: If „enable raster on kontur“ is selcted, this entry is editable. If the ratio of total time for the fast raster to total time for the normal raster is lower than the defined threshold in this entry, the fast raster will be acquired. If not, then only the normal raster will be acquired. Is state in percent and from 0 to 100.

center detection settings:

- **enable center detection during raster**: Enables the analysis of the center detection method from **Tools** during the acquisition of a raster, whereby the collected data is available direct after the finish of the acquisition. To guarantee that, the exposure time of one image should be greater than the process time of the method for one picture.
- **search width**: Defines the width of the search area in pixel. Has to be even to guarantee a symmetric dispersion relative to the center of the search area. $\text{width} \in \{1,975\} \subset \mathbb{N}$

- **search height**: Defines the height of the search area in pixel. Has to be even to guarantee a symmetric dispersion relative to the center of the search area. $\text{height} \in \{1,643\} \subset \mathbb{N}$
- **search x-offset**: Defines the offset in x-direction (cartesian coordinate system) in pixel from the center of the image to the center of the search area. $\text{x-offset} \in \{1,487\} \subset \mathbb{N}$
- **search y-offset**: Defines the offset in y-direction (cartesian coordinate system) in pixel from the center of the image to the center of the search area. $\text{y-offset} \in \{1,321\} \subset \mathbb{N}$
- **search linewidth**: Defines the half of the linewidth of each line in a mask in pixel. Values greater than 100 should be avoided, because otherwise a too small area is left to calculate the results and the results aren't significant. $\text{linewidth} \in \mathbb{N}$
- **starting angle**: Defines by which angle the first mask is rotated before the other masks are generated. It's important for the case that one „line“ of the diffraction pattern is nearly vertical, so a good boundary isn't visible any more in the visualisation of the rotated angle. $\text{angle} \in \{0,360\} \subset \mathbb{N}$
- **number of lines**: Defines the number of lines in one mask, which are all rotated by equally spaced angles. Therefor only numbers up to 10 are supported. $\in \{1,10\} \subset \mathbb{N}$
- **save info image during center detection**: Defines if an info-image is saved additional to the other acquired data. The info-image contains a visualisation of the search field and the detected center on the image and the information of the rotated angle and coordinates from the center.

apply:

Saves the changes and inspect all values for format errors. If a format or value error is arised, the settings dialog will not be closed and the corresponding error is print at [Frame 3](#)

default:

All entries are set to their default values. To change the default values, the user can modify the corresponding textfile (see also fig. 1.8) in the application folder where all default values are saved. The textfile is formatted in lines according to the order in the settings frame. Name and value of the parameter have to be separated by „:“ (without space!) and the value has to be in the same format like in the settings. The values of the checkboxes are 0 for not ticked and 1 for ticked.

```
1 gamma max:2.0
2 step size max:10.0
3 angle corr max:2.0
4 record 12-Bit images:0
5 scroll faktor:1.0
6 scroll step size int:1
7 scroll step size float:0.01
8 symmetry center x:487.5
9 symmetry center y:321.5
10 sample name:Sample_1
11 voltage:20
12 current:30
13 distance:40.0
14 comment:comment
15 normal bg-remove:0
16 remove dark-current:1
17 auto delta for polyfit:40
18 auto increment for polyfit:40
19 auto bg-remove during raster:0
20 calculate max exposure during raster:1
21 enable raster on kontur:0
22 threshold for kontur-raster:5.0
23 center detection during raster:0
24 search width for center detection:41
25 search height for center detection:41
26 search x offset for center detection:0
27 search y offset for center detection:0
28 linewidth for center detection:10
29 starting angle for center detection:0
30 number of lines:6
31 save info images during center detection:0
```

Fig. 1.8: default settings of the default settings in the textfile.

Frame 2 - image viewer

Frame 2 is a 8-bit grayscale image viewer, who can display images with 975x643 pixel. All modification made in [Section 4](#) will be immediately visible in this frame.

Frame 3 - console

In frame 3, important messages for the user will be displayed, so no editing is possible. If the console doesn't show any text, a simple mouseclick onto the frame or scroll the mousewheel is normally enough to update the console, if the mouse is in the frame.

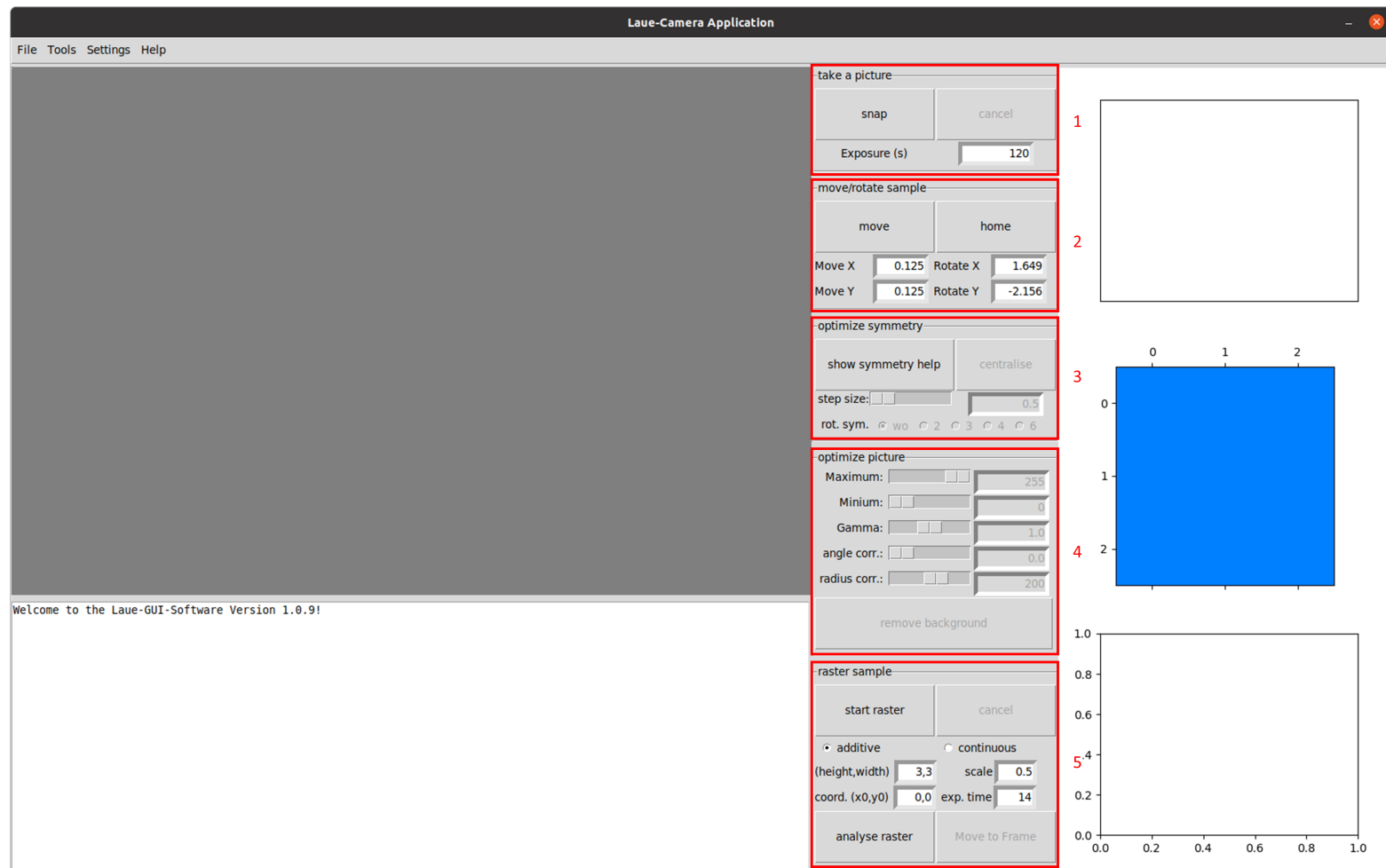


Fig. 4.1: Frame 4 divided into sections

Frame 4 - area of operations

Section 1 - image acquisition

snap:

An image with the defined exposure time will be acquired and afterwards displayed at [Frame 2](#), also the values of [Section 4](#) will be updated. Unfortunately the detector contains only a built-in 12-bit grayscale converter, so the resulting image is limited to 12-bit (0-4095). Therefore for all images with exposure times greater than 90 s (empirical value), the specific maximum exposure time has to be calculated. To achieve this, a separate image with 30 s exposure time will be acquired. With the maximum pixel value g_{30s} of this image and due to a linear correlation of maximum pixel value g_{\max} and exposure time, the maximum exposure time t_{\max} can be approximated with:

$$t_{\max} = \frac{g_{\max} - 100}{g_{30s} - 100} \cdot 30 \text{ s.} \quad (4.1)$$

The necessary offset of 100 results from the mean background noise. The maximum pixel value is set to 3900, so that the approximation doesn't reach by accident the limit value. If the calculation is completed, the defined exposure time will be divided into frames with the maximum exposure time plus one frame with the remaining exposure time. Afterwards all frames will be added to form the actual image.

cancel:

Cancels the actual image acquisition and reset the displayed image at [Frame 2](#).

Exposure (s):

Defines the exposure time. $\text{exp} \in \mathbb{N}$

Section 2 - translation and rotation of the sample

move:

The goniometer will be moved and rotated to the defined absolute values. To avoid the critical areas of the goniometer, only translations to ± 15 mm and rotations to $\pm 19^\circ$ are possible. To ensure this, it is recommend to manually reorientate the goniometer every time at the beginning of a measurement with a

new sample. Specifically the user has to push the **home**-Button, then manually adjust the rotators into their null positions and adjust the sample, so that the x-ray beam hits the center of the sample. Also it is important, that the coordinate system of the goniometer isn't a normal Cartesian coordinate system (see fig. 4.1).

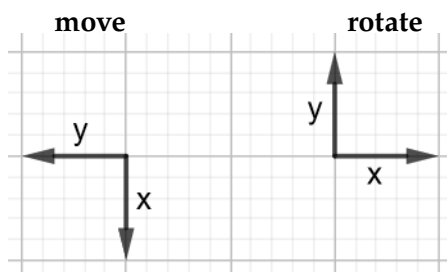


Fig. 4.1: coordinate systems of the goniometer.

home:

Reset all values to 0 and moves/rotates the goniometer in null position.

Move X:

Moves to the absolute value in x direction of the goniometer coordinate system, is state in mm. $\text{move}_x \in (-15.0, 15.0)$

Move Y:

Moves to the absolute value in y direction of the goniometer coordinate system, is state in mm. $\text{move}_y \in (-15.0, 15.0)$

Rotate X:

Rotates to the absolute value in x direction of the goniometer coordinate system, is state in °. $\text{rotate}_x \in (-19.0, 19.0)$

Rotate Y:

Rotates to the absolute value in y direction of the goniometer coordinate system, is state in °. $\text{rotate}_y \in (-19.0, 19.0)$

Section 3 - symmetrization

show symmetry help:

A red dot appears in [Frame 2](#) at the defined symmetry position in [Settings](#). If in „**rot. sym.**“ a symmetry is additionally selected, a vertical line from the dot to the upper edge is drawn plus the corresponding lines of the selected symmetry. The dot is movable with the arrow keys in „**step size**“-pixel steps and if the user clicks with the left mouse button in [Frame 2](#), the dot appears at this position. If the symmetry help lines are displayed, the user can rotate each

of these lines by „step size“-degree by scrolling with the mousewheel. This method offers the user on the one hand the opportunity to measure the angle between two points respectively to another point and on the other hand to centralize a reflex in the dead zone of the detector (see also „centralise“). If symmetry help is active, the user can disable this method by pushing the button again.

centralise:

This method offers the user to centralize the selected point via „show symmetry help“ on the symmetry coordinates. The main reason for this is to hide a bright reflex at the dead area of the detector, so minor reflexes are easier to detect. If the certain reflex is already inside the dead area, the lines from symmetry help allows a simpler way of detection.

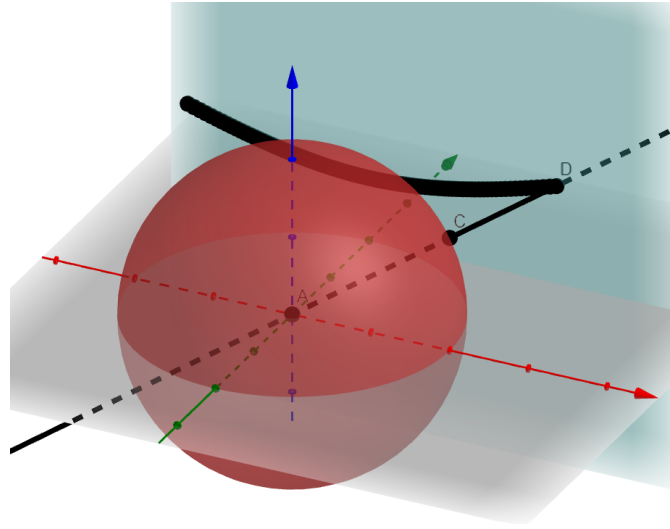


Fig. 4.2: Model to calculate the necessary rotational angles to centralize a reflex at a certain point at the detector. Moreover the variation of ϕ is displayed, to visualize the ϕ -dependency of θ .

The model for the necessary calculation is shown in fig. 4.2 with the target in the origin of the coordinate system and the detector at a distance d in xz -plane. With the help of a spherical coordinate system, the point of intersection D of line AC with the detector can be expressed and thus also the angles θ und ϕ by:

$$\phi = \arctan\left(\frac{d}{x}\right) \quad (4.2)$$

$$\theta = \arctan\left(\frac{d}{y \cdot \sin(\phi)}\right) \quad (4.3)$$

With both angles, the goniometer has to be rotated by these negated values to centralize the defined point to the center. Adapted to the coordinate system of

the goniometer and respectively to the current rotational values r_x and r_y , the new values can be determined with:

$$\text{rotate}_x = r_x - \frac{(90^\circ - \phi)}{2} \quad (4.4)$$

$$\text{rotate}_y = r_y - \frac{\theta - (90^\circ)}{2} \quad (4.5)$$

step size:

Defines the step size for the translation of the red dot of „show symmetry help“ in pixel and the rotation of the symmetry help lines in degree by scrolling the mousewheel. Is always rounded to the next number divisible by 0.5. $\text{stepsize} \in \mathbb{R}^+$

rot. sym.:

Defines the number of symmetry help lines. Selectable are values between 2 and 6 or none (wo).

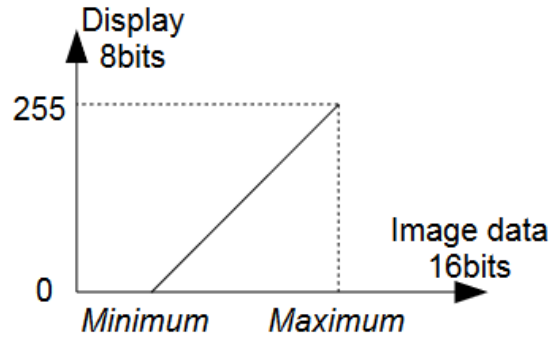


Fig. 4.3: Transformation between the selected range of minimum and maximum using the example of a transformation from 16-bit to 8-bit.

Section 4 - optimize picture

Maximum:

Defines the maximal value of the grayscale image, which is displayed in [Frame 2](#). All greater values are set to this maximum value. The selected interval between minimum and maximum is linear converted to 8-bit to display the resulting image at [Frame 2](#) (see [fig. 4.3](#)). Selectable range is between minimum and maximum pixel value. $\text{max} \in \mathbb{N}$

Minimum:

Defines the minimal value of the grayscale image, which is displayed in [Frame 2](#). All smaller values are set to this minimum value. The selected interval between minimum and maximum is linear converted to 8-bit to display the resulting image at [Frame 2](#) (see [fig. 4.3](#)). Selectable range is between minimum pixel value and maximum. $\min \in \mathbb{N}$

gamma:

Defines the gamma exponent, if a non linear contrast optimization is desired. Is applied on the image A after the definition of minimum and maximum, but before the transformation to an 8-bit image B , which is displayed at [Frame 2](#). The mathematical expression is:

$$B = A^\gamma \cdot 255^{1-\gamma} \quad (4.6)$$

angle corr.:

!!!Attention!!!: this and the following method are not longer necessary, but are not deleted if another purpose will show up. This optimization is only available after the background subtraction to reduce the remaining background noise. This function is only a experimentally method, which changes the data, so the usage must be marked and restricted to visual optimization, it is not useful for intensity comparison.

The remaining background noise is only located in the border area of the picture, so it will be concentrated on this area. The main idea is, that all reflexes will appear at the same point, but the background noise will be randomly spread. If the diffraction image is rotated around the center by small angles and subtracted from the original image, the noise will average out partially but the subtraction of the noise from the reflex will not affect it very much. Also several negative values will appear, which can be set to 0. Furthermore the image has to be rotated in both directions and weight half to avoid artefacts. This method will be applied to the image after minimum, maximum and gamma parameters are evaluated and the image is an 8-bit image. For a better visualization, the mean value of the picture will be subtracted and all negative values are set to 0 plus a transformation to the complete range of an 8-bit image. Selectable are rotation angles between 0 and the defined maximum value in [Settings](#).

radius corr.:

Because the method [angle corr.](#) is only suitable for border areas and arise false maxima in the center, this method enables to define the radius of a circle around the center in pixel, which displays the original image. The resulting image is way better than the former one, although the transition between the original image and the calculated image is visible. Selectable are radii between 0 and 321 (half height of the image). $\text{radius corr.} \in \mathbb{N}^+$

remove background:

Removes the background of the image. Therefore a central pixel line is selected from the original image, to which data the function

$$I = I_0 \cdot \left[\arctan \left(\frac{d - a \cdot (x - 487.5)}{g} \right) \arctan \left(\frac{-d - a \cdot (x - 487.5)}{g} \right) \right] + c \quad (4.7)$$

for the pixel values (I) depending on the distance in pixels from the left edge of the image (x) is fitted. This fitfunction is from a model, where an emitter with a $2d$ spread (depending on the aperture) is at the distance g to the detector. The parameter a represent the ratio between distance and pixel density, thus the horizontal spread of one pixel. The parameter d is set so $d = 40.2$, otherwise the fitfunction converges not every time and if it converges, no noteworthy differences are visible. Therefore it is recommend to choose approximately 40 mm as the distance between sample and detector. Other models, for examples the detector is angled relative to the x-ray beam or a round surface-emitter, are taken into account, but the resulting fitfunctions only converge very infrequently and if it converges, again no noteworthy differences are visible. On the other hand the asymmetry of both detector halves can be corrected easily. Therefore again a centered pixel line is selected, but this time the line is divided in the middle and both sectors are fitted separately with the fitfunction. Also both resulting fitfunction are subtracted separately from the image. The differences of both fitmodi can be seen in fig. 4.4.

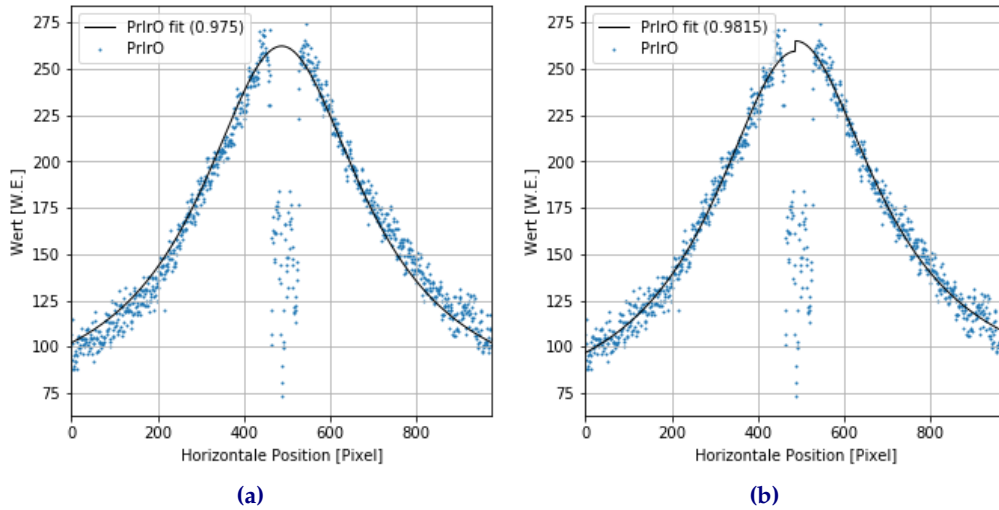


Fig. 4.4: Different results of both fitmodi. In (a) a normal fit with the data of the whole line is displayed and in (b) a partial fit for every detector half separately is displayed.

The passed image to the function is the original picture after the evaluation of

minimum, maximum and gamma exponent, but before the transformation to an 8-bit grayscale image. Therefore the return image is also a grayscale image with greater values than 8-bit. The original image with removed background can also be saved by [Save Original](#). After the background subtraction all image processing methods are available again plus „[angle corr.](#)“ and „[radius corr.](#)“ (**Attention:** usually not longer necessary).

Despite the background remove another background remains (see fig. 4.6c), which can be lead back to a dark-current inside the detector. To minimize this background, another method must be used, in which also the whole image is passed to the correspondig method if it is in [Settings](#) ticked. First, all pixel values greater than a maximum value are sorted out, so that the profile only contains the background and no great maxima. This maximum value is either calculated or determined manually with the [bg fine-tuning](#)-button. For the automatic calculation, the image is divided into the two detector halves and the mean value of every line of both sections is calculated, but at both sides of every line δ -pixel are excluded to avoid bad results due to the borders and the dead area in the center. The maximum value is set to the greates mean value plus an increment, which can be defined in [Settings](#) with δ . In fig. 4.5a a resulting line profile is shown exemplary. Because the profile is fluctuating very much, a 3x3-median filter is applied, so that the orange dots result in fig. 4.5b. Finally for each line in each detector half a polynomial function of 10th degree (blue line) is fitted and subtracted from the corresponding line. At the end in fig. 4.6d the result of this correction can be seen.


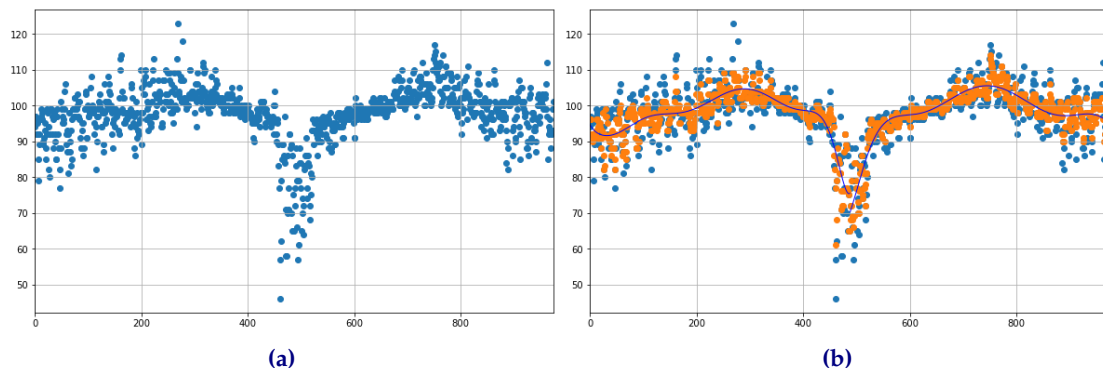


Fig. 4.5: In (a) a resulting line profile from the dark-current background is shown exemplary. In (b) a 3x3-median filter is applied onto the whole image and for each half a polynomial of 10th degree is fitted. This fitfunction describes the profile relatively good and is below the weak maxima, so that these maxima won't be immediately subtracted with the noise.


bg fine-tuning:

Because the automatic calculation of the maximum pixel value is good but not perfect all the time, this method gives the user a opportunity to adjust the


value manually. In the opening window the actual image without background is displayed, but all pixel values greater than the defined maximum value at the slider aren't displayed. With the slider, the user can vary the maximum value and see the resulting image immediately. In general only the reflexes should be excluded and the background from the dark-current should remain. With the **apply**-button, the new maximum value is set to the selected value and the displayed image in the main window is updated. With the **default**-button the slider is set back to the calculated value.

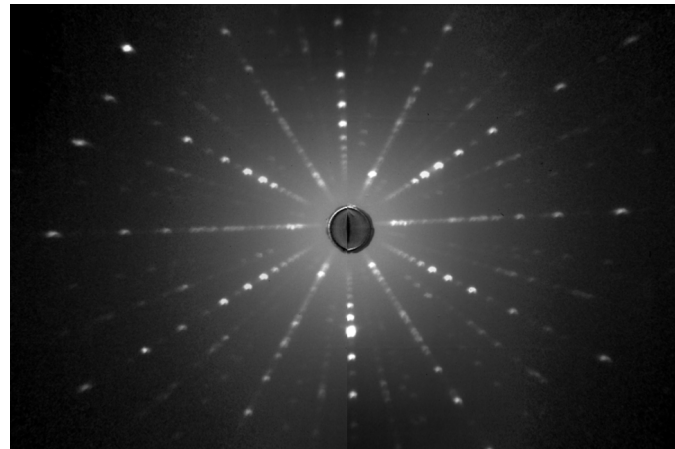


(a) original image

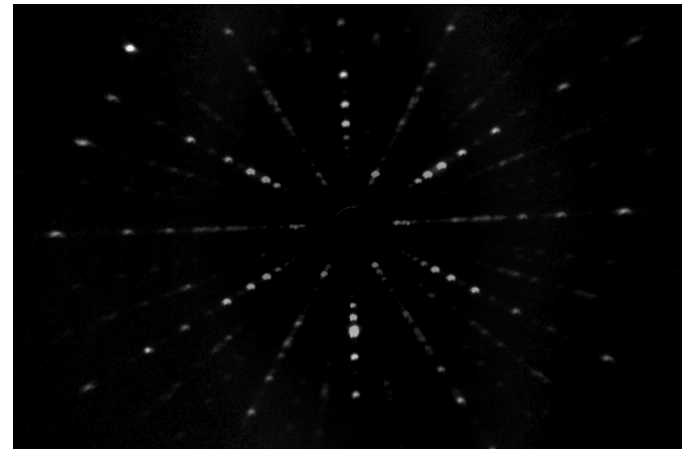

(b) image without background

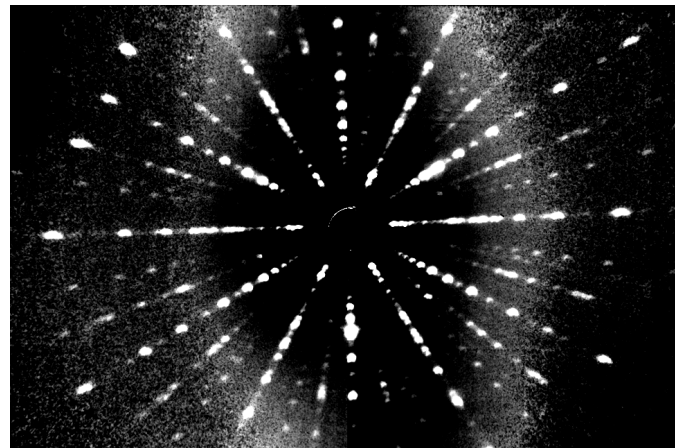

(c) image without background and with contrast optimization

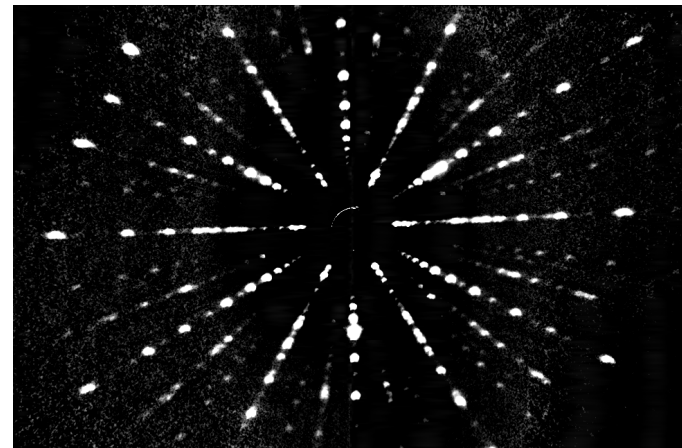

(d) image without background, with contrast optimization plus polynomial fitting and so without dark-current background

Fig. 4.6: Demonstration of a possible image processing of a original laue-image

Section 5 - raster sample

start raster:

Starts a new raster with the defined values and the selected settings in [Settings](#). Furthermore for each acquired frame a classification is made on how good the diffraction pattern is (from 0 (not visible) to 2 (good)). In the contour diagram in [Frame 5](#) the actual progress is displayed and in [Frame 2](#) the last acquired image is visible. At this image also the methods from [Section 4](#) are applicable. In the beginning of the acquisition of a raster, a new folder is created at the selected directory with the filename „name/number of sample“ - „voltage“ - „current“ - „distance to detector in mm“ - „RASTER“ - „exposure time“ - „height x width“ - „stepsize in mm“ - „additive or continuous“. Additionally at the end of a raster a „logfile.txt“ is created which contains all information of the raster settings and the acquired contour (for the format see a corresponding logfile.txt). With this logfile also raster which were acquired before this application was written can be displayed by creating a logfile.txt by hand and use „change kontur manually“ to change the contour.

cancel:

Cancels the actual raster.

additive, continuous:

Defines if a frame is acquired with a maximum exposure time of 16 s (additive) or with the calculated maximum exposure time (continuous). The advantage of additive over continuous is a better contour detection of the algorithm, but the advantage of the continuous method is a longer uninterrupted exposure time, thus weak reflexes can be seen better.

height, width:

Defines the height and width of the raster, thus how many single frames are acquired by each row and column. $h, w \in \mathbb{N}^+, \mathbb{N}^+$

scale:

Defines the step size in mm between two frames in a raster. $scale \in \mathbb{R}^+$

coord. (x_0, y_0) :

Defines the center of the raster, from which the corresponding raster will start. $x_0, y_0 \in \mathbb{R}^+, \mathbb{R}^+$

exp. time:

Defines the exposure time in seconds of each frame. $exposure \in \mathbb{R}^+$

analyse raster:

If a raster is acquired or loaded, with this button the particular contour can be analysed. It gives the user the opportunity to navigate with the arrow keys

through the raster. Furthermore the current frame is visualized through a brown frame and the corresponding diffraction image is shown in [Frame 2](#), ready to be modified by the optimization methods from [Section 4](#). Immediately after a raster acquisition, the contour is deposit in [Frame 5](#).

The contour navigation can be canceled by clicking the button again.

move to frame:

If the current raster and the current coordinate system of the goniometer are the same, the position of the sample can be set to the selected frame of the contour.

Frame 5 - diagram viewer

Diagram 1 - maximum finder

In this diagram all positions of the detected maxima of the current frame in an acquired raster is displayed. With this additional information the user can confirm or falsify the classification done by the algorithm. If a multiplot raster is loaded, only the corresponding part of the image is displayed.

If an original raster is loaded, the maxima will be detected by the same function as in the multiplot but this method will first be applied if the current position of the contour changes.

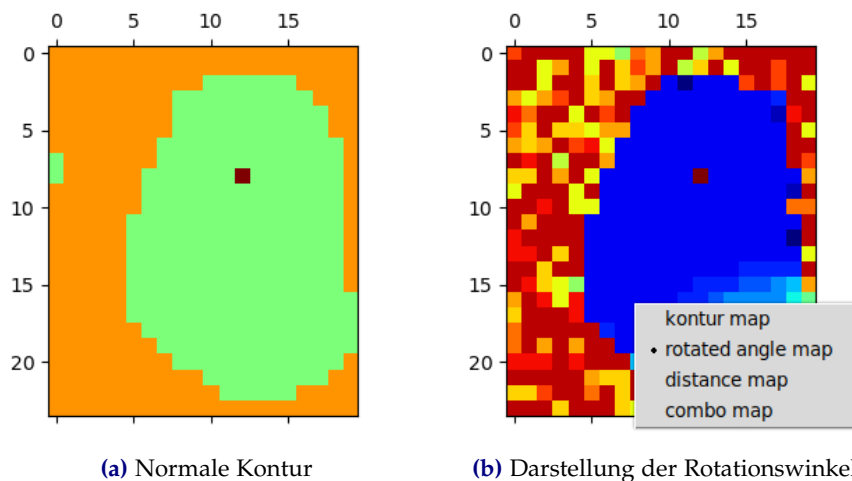


Fig. 5.1: Different types of the visualization of the kontur from a raster with center detection data.

Diagram 2 - contour diagram

In this diagram the current contour is displayed, which will show the current progress of a raster acquisition or the whole contour of an already measured raster. If the navigation on the raster is active, the current position is marked by a brown frame. If a raster with center detection data is loaded, these data can be visualized in this diagram with different types, selectable through a right-click on the diagram.

In the method „**auto find kontur**“, the contour will build up time after time, so the user can visually comprehend the contour detection.

Diagram 3 - backgroundfit

If the background is removed from an image, this diagram will visualize the fit with the data from the image (blue dots) and the fitfunction (black line). The left axis shows the pixel value of the corresponding pixel and the lower axis shows the horizontal pixel number from the left image border.

Outlook

Upcoming improvements:

- Another visualization of the contours, which takes into account that the detected area of the sample is circular and not quadratic and frames are overlapping with a small stepsize.
- Make possible to detect smaller differences in the diffraction pattern in a raster and display them.
- Drop out of a raster before it is completely acquired, if the sample is fully detected but several frames are in the line. The remaining frames will be filled up with standard images, so that an analysis with the application still will be possible.
- Automatising of the centralise-method for good diffraction images.
- Introduce another nonlinear optimization function for the contrast besides the gamma factor, to make weak reflexes better visible.

If unknown errors or other questions arise or more extensions are reasonable, please contact the author!